

EDELVITAL™

3 IN 1 RAPID ANTIGEN CORONA TEST

99,52%*

More on:
ameyamed.com



MULTI PACK



EDELVITAL™ 3 IN 1 RAPID ANTIGEN TEST KIT (NASAL, THROAT, SALIVA)

		Antigen test 3 in 1
Test type		Front nasal, throat, saliva
Accuracy		99,52%*
Sensitivity		99,18%*
Specificity		100%*
Durability		24 Months
Storage temp.		2-30°C
PU		20 Pcs. per kit

The clinical performance of the rapid COVID-19 antigen test (colloidal gold) was determined by testing 1096 positive and 793 negative samples for SARS-CoV-2 antigen with a sensitivity of 99.18% (95% CI: 98.45% - 99.62%) and a specificity of 100% (95% CI: 99.54% - 100%). Clinical specimens were determined to be positive or negative using an RT-PCR reference method.

THE 3IN1 RAPID ANTIGEN TEST FOR EASIEST HANDLING IN JUST 3 WORK STEPS

The EDELVITAL™ Rapid COVID-19 antigen test is a colloidal gold immunochromatography for the qualitative detection of nucleocapsid antigens from SARS-CoV-2 in human nasal swabs, throat swabs or saliva from people suspected of having COVID-19.

The results relate to the identification of the SARS-CoV-2 nucleocapsid antigen. The antigen is generally detectable in samples from the upper or lower respiratory tract during the acute phase of infection. The positive results indicate the presence of viral antigens, but clinical correlation with medical history and other diagnostic information is required to determine infection status.

THE BENEFITS



- 3 in 1 specimen collection with front nasal, throat or comfortable saliva swab.
- High performance proven by clinical study.
- Fast detection of SARS-COV-2 nucleocapsid antigen in 15 minutes.
- No instrumentation required.
- Inexpensive solution for large-scale testing.

- 2-30° Storage
- 15 MIN Duration
- 99,18%* Sensitivity
- 100%* Specificity
- 99,52%* Accuracy



**INQUIRY
PACKAGING UNITS
UNIT PRICES & PALLET GOODS**

Packaging unit:
20 pieces per pack,
500 pieces per box

Priority delivery and business inquiry please via
phone: **+43 1 523 3000** or email:
info@ameyamed.com

Ameya Med is a brand of the Ameya GmbH.
Address: Ameya GmbH, Gregor-Mendel-Strasse
46, 1190 Vienna, Austria

Tel.: +43 1 523 3000
Fax: +43 1 523 3000 33

Website: www.ameyamed.com
Email: info@ameyamed.com


VAT: ATU63103714
CRN: 285467f
GISA No.: 33180168:

**1 PKG TEST KIT
CONTAINS:**

- 20x Test card
- 20x Extractiontube with pre-filled solution
- 20x Swabs
- 1x User manual

**TEST
PROCEDURE**

Read the instructions for use
for the exact process.

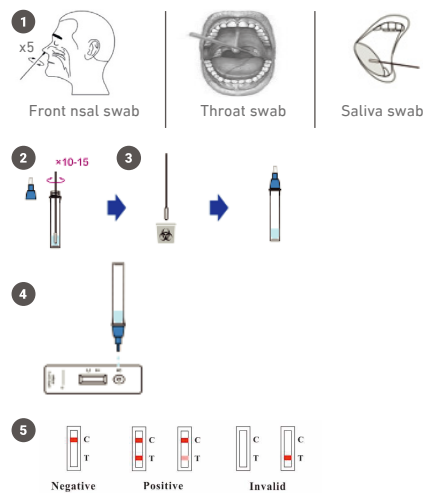

HIGHEST
QUALITY